

# SVENSK STANDARD

## SS-EN 15703-2:2014

Fastställt/Approved: 2014-12-21  
Publicerad/Published: 2015-01-09  
Utgåva/Edition: 1  
Språk/Language: engelska/English  
ICS: 77.040.30; 77.120.30

---

### **Koppar och kopparlegeringar – Bestämning av manganhalt – Del 2: Flamatomär absorptionspektrometrimetod (FAAS)**

**Copper and copper alloys – Determination of manganese  
content –  
Part 2: Flame atomic absorption spectrometric method (FAAS)**

This preview is downloaded from [www.sis.se](http://www.sis.se). Buy the entire standard via <https://www.sis.se/std-104944>

# Standarder får världen att fungera

*SIS (Swedish Standards Institute) är en fristående ideell förening med medlemmar från både privat och offentlig sektor. Vi är en del av det europeiska och globala nätverk som utarbetar internationella standarder. Standarder är dokumenterad kunskap utvecklad av framstående aktörer inom industri, näringsliv och samhälle och befrämjar handel över gränser, bidrar till att processer och produkter blir säkrare samt effektiviserar din verksamhet.*

## Delta och påverka

Som medlem i SIS har du möjlighet att påverka framtida standarder inom ditt område på nationell, europeisk och global nivå. Du får samtidigt tillgång till tidig information om utvecklingen inom din bransch.


## Ta del av det färdiga arbetet

Vi erbjuder våra kunder allt som rör standarder och deras tillämpning. Hos oss kan du köpa alla publikationer du behöver – allt från enskilda standarder, tekniska rapporter och standardpaket till handböcker och onlinetjänster. Genom vår webbtjänst e-nav får du tillgång till ett lättnavigerat bibliotek där alla standarder som är aktuella för ditt företag finns tillgängliga. Standarder och handböcker är källor till kunskap. Vi säljer dem.

## Utveckla din kompetens och lyckas bättre i ditt arbete

Hos SIS kan du gå öppna eller företagsinterna utbildningar kring innehåll och tillämpning av standarder. Genom vår närhet till den internationella utvecklingen och ISO får du rätt kunskap i rätt tid, direkt från källan. Med vår kunskap om standarders möjligheter hjälper vi våra kunder att skapa verklig nytta och lönsamhet i sina verksamheter.

**Vill du veta mer om SIS eller hur standarder kan effektivisera din verksamhet är du välkommen in på [www.sis.se](http://www.sis.se) eller ta kontakt med oss på tel 08-555 523 00.**

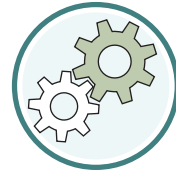

# Standards make the world go round

*SIS (Swedish Standards Institute) is an independent non-profit organisation with members from both the private and public sectors. We are part of the European and global network that draws up international standards. Standards consist of documented knowledge developed by prominent actors within the industry, business world and society. They promote cross-border trade, they help to make processes and products safer and they streamline your organisation.*

## Take part and have influence

As a member of SIS you will have the possibility to participate in standardization activities on national, European and global level. The membership in SIS will give you the opportunity to influence future standards and gain access to early stage information about developments within your field.

## Get to know the finished work

We offer our customers everything in connection with standards and their application. You can purchase all the publications you need from us - everything from individual standards, technical reports and standard packages through to manuals and online services. Our web service e-nav gives you access to an easy-to-navigate library where all standards that are relevant to your company are available. Standards and manuals are sources of knowledge. We sell them.

## Increase understanding and improve perception

With SIS you can undergo either shared or in-house training in the content and application of standards. Thanks to our proximity to international development and ISO you receive the right knowledge at the right time, direct from the source. With our knowledge about the potential of standards, we assist our customers in creating tangible benefit and profitability in their organisations.

**If you want to know more about SIS, or how standards can streamline your organisation, please visit [www.sis.se](http://www.sis.se) or contact us on phone +46 (0)8-555 523 00**



Europastandarden EN 15703-2:2014 gäller som svensk standard. Detta dokument innehåller den officiella engelska versionen av EN 15703-2:2014.

The European Standard EN 15703-2:2014 has the status of a Swedish Standard. This document contains the official version of EN 15703-2:2014.

© Copyright/Upphovsrätten till denna produkt tillhör SIS, Swedish Standards Institute, Stockholm, Sverige. Användningen av denna produkt regleras av slutanvändarlicensen som återfinns i denna produkt, se standardens sista sidor.

© Copyright SIS, Swedish Standards Institute, Stockholm, Sweden. All rights reserved. The use of this product is governed by the end-user licence for this product. You will find the licence in the end of this document.

*Uppllysningar om sakinnehållet i standarden lämnas av SIS, Swedish Standards Institute, telefon 08-555 520 00. Standarder kan beställas hos SIS Förlag AB som även lämnar allmänna uppllysningar om svensk och utländsk standard.*

*Information about the content of the standard is available from the Swedish Standards Institute (SIS), telephone +46 8 555 520 00. Standards may be ordered from SIS Förlag AB, who can also provide general information about Swedish and foreign standards.*

Denna standard är framtagen av kommittén för Koppar, SIS/TK 132.

Har du synpunkter på innehållet i den här standarden, vill du delta i ett kommande revideringsarbete eller vara med och ta fram andra standarder inom området? Gå in på [www.sis.se](http://www.sis.se) - där hittar du mer information.



EUROPEAN STANDARD

**EN 15703-2**

NORME EUROPÉENNE

EUROPÄISCHE NORM

December 2014

---

ICS 77.040.30; 77.120.30

English Version

## Copper and copper alloys - Determination of manganese content - Part 2: Flame atomic absorption spectrometric method (FAAS)

Cuivre et alliages de cuivre - Détermination de manganèse  
- Partie 2: Méthode par spectrométrie d'absorption  
atomique dans la flamme (SAAF)

Kupfer und Kupferlegierungen - Bestimmung des  
Mangangehaltes - Teil 2:  
Flammenatomabsorptionsspektrometrisches Verfahren  
(FAAS)

This European Standard was approved by CEN on 8 November 2014.

CEN members are bound to comply with the CEN/CENELEC Internal Regulations which stipulate the conditions for giving this European Standard the status of a national standard without any alteration. Up-to-date lists and bibliographical references concerning such national standards may be obtained on application to the CEN-CENELEC Management Centre or to any CEN member.

This European Standard exists in three official versions (English, French, German). A version in any other language made by translation under the responsibility of a CEN member into its own language and notified to the CEN-CENELEC Management Centre has the same status as the official versions.

CEN members are the national standards bodies of Austria, Belgium, Bulgaria, Croatia, Cyprus, Czech Republic, Denmark, Estonia, Finland, Former Yugoslav Republic of Macedonia, France, Germany, Greece, Hungary, Iceland, Ireland, Italy, Latvia, Lithuania, Luxembourg, Malta, Netherlands, Norway, Poland, Portugal, Romania, Slovakia, Slovenia, Spain, Sweden, Switzerland, Turkey and United Kingdom.



EUROPEAN COMMITTEE FOR STANDARDIZATION  
COMITÉ EUROPÉEN DE NORMALISATION  
EUROPÄISCHES KOMITEE FÜR NORMUNG

**CEN-CENELEC Management Centre: Avenue Marnix 17, B-1000 Brussels**

| <b>Contents</b>   | <b>Page</b> |
|---|-------------|
| Foreword.....   | 3           |
| 1 Scope .....   | 4           |
| 2 Normative references .....  | 4           |
| 3 Principle.....  | 4           |
| 4 Reagents.....   | 4           |
| 4.1 Hydrochloric acid, HCl ( $\rho = 1,19$ g/ml).....                         | 4           |
| 4.2 Nitric acid, HNO <sub>3</sub> ( $\rho = 1,40$ g/ml).....                  | 4           |
| 4.3 Hydrofluoric acid, HF ( $\rho = 1,13$ g/ml).....                          | 4           |
| 4.4 Nitric acid solution, 1 + 1 .....   | 4           |
| 4.5 Lanthanum(III) chloride solution, 100 g/l .....                           | 4           |
| 4.6 Sulphuric acid, H <sub>2</sub> SO <sub>4</sub> ( $\rho = 1,84$ g/ml)..... | 4           |
| 4.7 Sulphuric acid solution, 1 + 9.....                                       | 5           |
| 4.8 Manganese stock solution, 1,0 g/l Mn .....                                | 5           |
| 4.9 Manganese standard solution, 0,10 g/l.....                                | 5           |
| 4.10 Manganese standard solution, 0,010 g/l.....                              | 5           |
| 4.11 Copper matrix solution, 20 g/l Cu .....                                  | 5           |
| 4.12 Copper matrix solution, 2 g/l Cu .....                                   | 5           |
| 5 Apparatus .....   | 5           |
| 5.1 Atomic absorption spectrometer, fitted with an air/acetylene burner ..... | 5           |
| 5.2 Manganese hollow-cathode lamp .....                                       | 6           |
| 6 Sampling.....   | 6           |
| 7 Procedure .....   | 6           |
| 7.1 Preparation of the test portion solution .....                            | 6           |
| 7.2 Blank test.....   | 7           |
| 7.3 Check test.....   | 7           |
| 7.4 Establishment of the calibration curve .....                              | 7           |
| 7.5 Determination.....  | 10          |
| 8 Expression of results .....   | 10          |
| 8.1 Use of the calibration curve.....   | 10          |
| 8.2 Use of the bracketing method .....  | 11          |
| 9 Precision.....  | 12          |
| 10 Test report .....  | 13          |
| Bibliography .....  | 14          |

## Foreword

This document (EN 15703-2:2014) has been prepared by Technical Committee CEN/TC 133 “Copper and copper alloys”, the secretariat of which is held by DIN.

This European Standard shall be given the status of a national standard, either by publication of an identical text or by endorsement, at the latest by June 2015 and conflicting national standards shall be withdrawn at the latest by June 2015.

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. CEN [and/or CENELEC] shall not be held responsible for identifying any or all such patent rights.

Within its programme of work, Technical Committee CEN/TC 133 requested CEN/TC 133/WG 10 “Methods of analysis” to prepare the following standard:

EN 15703-2, *Copper and copper alloys — Determination of manganese content — Part 2: Flame atomic absorption spectrometric method (FAAS)*

This is one of two Parts of the standard/Technical Specification for the determination of manganese content in copper and copper alloys. The other Part is:

CEN/TS 15703-1, *Copper and copper alloys — Determination of manganese content — Part 1: Spectrophotometric method*

According to the CEN-CENELEC Internal Regulations, the national standards organizations of the following countries are bound to implement this European Standard: Austria, Belgium, Bulgaria, Croatia, Cyprus, Czech Republic, Denmark, Estonia, Finland, Former Yugoslav Republic of Macedonia, France, Germany, Greece, Hungary, Iceland, Ireland, Italy, Latvia, Lithuania, Luxembourg, Malta, Netherlands, Norway, Poland, Portugal, Romania, Slovakia, Slovenia, Spain, Sweden, Switzerland, Turkey and the United Kingdom.

## 1 Scope

This European Standard specifies a flame atomic absorption spectrometric method (FAAS) for the determination of the manganese content of copper and copper alloys in the form of unwrought, wrought and cast products.

The method is applicable to products having manganese mass fractions between 0,001 0 % and 6,0 %.

## 2 Normative references

The following documents, in whole or in part, are normatively referenced in this document and are indispensable for its application. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 1811-1, *Copper and copper alloys — Selection and preparation of samples for chemical analysis — Part 1: Sampling of cast unwrought products*

ISO 1811-2, *Copper and copper alloys — Selection and preparation of samples for chemical analysis — Part 2: Sampling of wrought products and castings*

## 3 Principle

Dissolution of a test portion in hydrochloric acid and nitric acid followed, after suitable dilution and the addition of lanthanum chloride to mask the effect of interfering ions, by aspiration into an air/acetylene flame of an atomic absorption spectrometer. Measurement of the absorption of the 279,5 nm or 403,1 nm line emitted by a manganese hollow-cathode lamp.

## 4 Reagents

During the analysis, use only reagents of recognized analytical grade and only distilled water or water of equivalent purity.

**4.1 Hydrochloric acid**, HCl ( $\rho = 1,19$  g/ml)

**4.2 Nitric acid**, HNO<sub>3</sub> ( $\rho = 1,40$  g/ml)

**4.3 Hydrofluoric acid**, HF ( $\rho = 1,13$  g/ml)

**WARNING — Hydrofluoric acid is extremely irritating and corrosive to skin and mucous membranes producing severe skin burns which are slow to heal. In the case of contact with skin, wash well with water, apply a topical gel containing 2,5 % (mass fraction) calcium gluconate, and seek immediate medical treatment.**

**4.4 Nitric acid solution**, 1 + 1

Add 500 ml of nitric acid (4.2) to 500 ml of water.

**4.5 Lanthanum(III) chloride solution**, 100 g/l

Dissolve 100 g of lanthanum(III) chloride (LaCl<sub>3</sub> · 7H<sub>2</sub>O) in a 600 ml beaker with water and transfer the solution into a 1 000 ml one-mark volumetric flask. Dilute to the mark with water and mix well.

**4.6 Sulphuric acid**, H<sub>2</sub>SO<sub>4</sub> ( $\rho = 1,84$  g/ml)



#### 4.7 Sulphuric acid solution, 1 + 9

Add 20 ml of sulphuric acid (4.6) in 180 ml water.

#### 4.8 Manganese stock solution, 1,0 g/l Mn

##### 4.8.1 Cleaning of manganese metal

Transfer several grams of electrolytic manganese (purity > 99,9 %) into a 250 ml beaker containing about 150 ml of sulphuric acid solution (4.7). Stir, then allow the manganese to settle for several minutes. Decant and discard the sulphuric acid solution. Rinse several times with water and finally with acetone. Dry the metal for about 2 min at 100°C and cool in a desiccator.

##### 4.8.2 Preparation of stock solution

Weigh,  $(1 \pm 0,001)$  g, of manganese cleaned as in 4.8.1 and transfer it into a 250 ml beaker. Add 20 ml of hydrochloric acid (4.1) and 20 ml of the nitric acid solution (4.4). Cover with a watch glass and, if necessary, heat gently to assist dissolution. When dissolution is complete, allow to cool and transfer the solution quantitatively into a 1 000 ml one-mark volumetric flask, dilute to the mark with water and mix well.

1 ml of this solution contains 1,0 mg of Mn.

#### 4.9 Manganese standard solution, 0,10 g/l

Transfer 20,0 ml of manganese stock solution (4.8) into a 200 ml one-mark volumetric flask. Dilute to the mark with water and mix well.

Prepare this solution immediately prior to use.

1 ml of this solution contains 0,10 mg of Mn.

#### 4.10 Manganese standard solution, 0,010 g/l

Transfer 5,0 ml of manganese stock solution (4.8) to a 500 ml one-mark volumetric flask. Dilute to the mark with water and mix well.

Prepare this solution immediately prior to use.

1 ml of this solution contains 0,010 mg of Mn.

#### 4.11 Copper matrix solution, 20 g/l Cu

Transfer  $(10 \pm 0,01)$  g of manganese-free copper ( $\text{Cu} \geq 99,95\%$ ) into a 600 ml beaker. Add 100 ml of hydrochloric acid (4.1) and, cautiously, 100 ml of the nitric acid solution (4.4). Cover with a watch glass and heat gently until the copper has been completely dissolved, then heat up to the boiling point until the nitrous fumes have been expelled. Allow to cool and transfer the solution quantitatively into a 500 ml one-mark volumetric flask, dilute to the mark with water and mix well.

#### 4.12 Copper matrix solution, 2 g/l Cu

Transfer 50,0 ml of the copper matrix solution (4.11) into a 500 ml one-mark volumetric flask. Add 90 ml of hydrochloric acid (4.1) and 90 ml of the nitric acid solution (4.4). Dilute to the mark with water and mix well.

## 5 Apparatus

### 5.1 Atomic absorption spectrometer, fitted with an air/acetylene burner

## 5.2 Manganese hollow-cathode lamp

## 6 Sampling

Sampling shall be carried out in accordance with ISO 1811-1 or ISO 1811-2, as appropriate.

Test samples shall be in the form of fine drillings, chips or millings with a maximum thickness of 0,5 mm.

## 7 Procedure

### 7.1 Preparation of the test portion solution

#### 7.1.1 Test portion

Weigh ( $1 \pm 0,001$ ) g of the test sample.

#### 7.1.2 Test portion solution

Transfer the test portion (7.1.1) into a 250 ml beaker. Add 10 ml of hydrochloric acid (4.1) and 10 ml of the nitric acid solution (4.4). Cover with a watch glass and heat gently until the test portion is completely dissolved. Allow to cool. If undissolved matter remains, indicating the presence of silicon, filter the solution. Place the filter and retained residues in a platinum crucible and ash, taking care that the filter does not flame. Calcine at about 550 °C. Allow to cool and add 5 ml of hydrofluoric acid (4.3) and five drops of nitric acid (4.2). Evaporate to dryness and calcine again for several minutes at about 700 °C to completely volatilize the silicon. Cool, and then dissolve the residue with the least possible volume of nitric acid solution (4.4). Filter, if necessary, and add this filtrate quantitatively to the original filtrate.

#### 7.1.3 Manganese mass fractions between 0,001 % and 0,020 %

Transfer the dissolved test portion or the combined filtrates quantitatively into a 100 ml one-mark volumetric flask. Add 10 ml of the lanthanum(III) chloride solution (4.5), dilute to the mark with water and mix well.

NOTE The validation exercise of this method showed that for manganese mass fractions between 0,015 % and 0,020 % results are better when the determination is carried out following 7.1.3.

#### 7.1.4 Manganese mass fractions between 0,020 % and 0,50 %

Transfer the dissolved test portion or the combined filtrates quantitatively into a 100 ml one-mark volumetric flask. Dilute to the mark with water and mix well. Transfer 5,0 ml of this solution into a 100 ml one-mark volumetric flask. Add 4,5 ml of hydrochloric acid (4.1), 4,5 ml of the nitric acid solution (4.4) and 10 ml of the lanthanum(III) chloride solution (4.5). Dilute to the mark with water and mix well.

#### 7.1.5 Manganese mass fractions between 0,50 % and 2,0 %

Transfer the dissolved test portion or the combined filtrates quantitatively into a 100 ml one-mark volumetric flask. Dilute to the mark with water and mix well. Transfer 10,0 ml of this solution into a 100 ml one-mark volumetric flask. Add 9 ml of hydrochloric acid (4.1), 9 ml of the nitric acid solution (4.4) and 10 ml of the lanthanum(III) chloride solution (4.5). Dilute to the mark with water and mix well.

NOTE The validation exercise of this method showed that for manganese mass fractions near 0,5 % results are better when the determination is carried out following 7.1.5.

#### 7.1.6 Manganese mass fractions between 2,0 % and 6,0 %

Transfer the dissolved test portion or the combined filtrates quantitatively into a 100 ml one-mark volumetric flask. Dilute to the mark with water and mix well. Transfer 2,5 ml of this solution into a 100 ml one-mark

volumetric flask. Add 2,3 ml of hydrochloric acid (4.1), 2,3 ml of the nitric acid solution (4.4) and 10 ml of the lanthanum(III) chloride solution (4.5). Dilute to the mark with water and mix well.

**NOTE** The validation exercise of this method showed that for manganese mass fractions near 2,0 % results are better when the determination is carried out following 7.1.6.

## **7.2 Blank test**

Carry out a blank test simultaneously with the determination, following the same procedure and using the same quantities of all reagents as used for the determination, but substituting pure copper for the test portion.

## **7.3 Check test**

Make a preliminary check of the apparatus by preparing a solution of a reference material or a synthetic sample containing a known amount of manganese and of composition similar to the material to be analysed. Carry out the procedure specified in 7.5.

## **7.4 Establishment of the calibration curve**

### **7.4.1 Preparation of the calibration solutions**

#### **7.4.1.1 General**

In all cases, copper, chloride and nitrate concentrations, and acidity in the calibration solutions shall be similar to those of the test portion solutions.

The presence of copper in the calibration solutions compensates for chemical interaction effects of copper in the test portion solution. Normally no similar additions are required to compensate for the effect of alloying elements. If an alloying element is present in the material to be analysed in mass fraction > 10 %, an appropriate mass of this element shall be added to the calibration solutions. The volumes of copper matrix solution added (4.11 or 4.12) have been calculated to compensate for chemical interaction effects of copper in test solutions of copper or high-copper alloys. Overcompensation may occur if the same volumes are added when the test samples are copper-based alloys where the percentage of copper is lower. In these cases the volumes of copper matrix solution shall be decreased to match the copper content of the test sample in solution.

The manganese concentration of the calibration solutions shall be adjusted to suit the sensitivity of the spectrometer used, so that the curve of absorbance as a function of concentration is a straight line.

#### **7.4.1.2 Manganese mass fractions between 0,001 0 % and 0,020 %**

Into each of a series of seven 100 ml one-mark volumetric flasks, introduce the volumes of manganese standard solution (4.10) and copper matrix solution (4.11) shown in Table 1. Introduce also 10 ml of lanthanum(III) chloride solution (4.5). Dilute to the mark with water and mix well.